

MONTAGE[®]

SETTABLE BONE PASTE



FLOWABLE

Contains Hydroxyapatite/ β -TCP and 5% by weight Vitamin E Acetate

MONTAGE *FLOWABLE* is a self-setting (hardening), resorbable bone paste that hardens in minutes following application. When mixed together (no hand-mixing required), the components of the MONTAGE *FLOWABLE* device form a cohesive and flowable bone paste.

PASTE CARTRIDGE



MIX TIP

Static mix tip automatically mixes the material as it is extruded.

Each mix tip features a radiopaque ring providing high-visibility under radiographic imaging.

Total mix tip length: ~103mm.

Three (3) mix tips are included in the packaging.



APPLICATOR

Ergonomic design maximizes comfort and control for precision delivery of material

PRODUCT DETAILS

- MONTAGE *FLOWABLE* is prefilled to deliver a minimum of 8cc of material
- MONTAGE *FLOWABLE* can be used up to three times throughout a procedure by changing the mix tip

PRODUCT CODE	PRODUCT DESCRIPTION
OS-MON-1501FL	MONTAGE <i>FLOWABLE</i> 8cc (15g)

PREPARATION & APPLICATION TECHNIQUE GUIDE

STEP 1 PREPARE THE BONE

- MONTAGE FLOWABLE adheres best to bone that has been debrided to remove all undesired tissue (loose bone, periosteum, soft tissue) and any other foreign materials.
- Before opening and applying MONTAGE FLOWABLE, be sure to irrigate and dry to ensure a clean bone surface.

STEP 2 ASSEMBLE THE DEVICE

- Open MONTAGE FLOWABLE packaging immediately prior to use.
- Attach the cartridge to the applicator ensuring that the larger chamber is oriented to mate with the larger of the two plungers on the applicator. **A B**

STEP 3 PRIME THE CARTRIDGE

The cartridge must be primed to balance the two paste components.

- Remove the cartridge cap by rotating it counter-clockwise and pulling it away from the cartridge.
- Balance the pistons in the cartridge by slowly squeezing the trigger of the applicator to express a small amount from both chambers. **C**
- Wipe any excess from the end of the cartridge.

STEP 4 APPLY



- Attach a static mix tip to the cartridge by aligning the triangular-tab on the mix tip with the corresponding notch on the cartridge. **D1**
- Press the mix tip firmly against the cartridge housing and rotate one quarter turn clockwise until it meets its locking stop. **D2**
- Slowly squeeze the trigger to express MONTAGE FLOWABLE from the end of the mix tip, to deliver the paste directly to the bone surface.
- MONTAGE FLOWABLE can be manipulated and spread over the bone surface following application with gloved fingers or a blunt surgical instrument with a working end that is larger than the intended area of application.

***Note: Squeezing harder will not extrude the paste quicker.**

Mix tips must be changed between applications as needed as MONTAGE FLOWABLE will harden inside the tip as soon as the dispensing process is halted.

STEP 5 REMOVE EXCESS

- Remove any excess MONTAGE FLOWABLE that is not necessary and/or impinging on adjacent soft tissues.
- Discard the device.

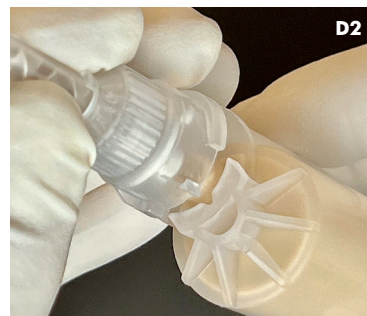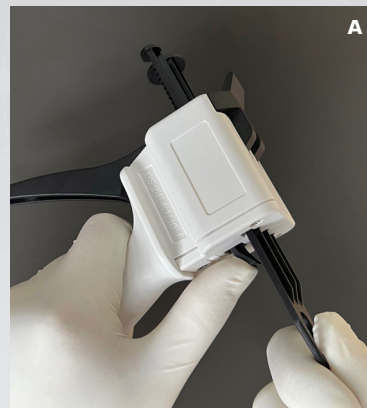

700 Fairfield Avenue, Suite 1 • Stamford, CT 06902 U.S.A. • 1.855.475.9175 • customerservice@abyrx.com • www.abyrx.com

© 2022 ABYRX, Inc.

ABYRX products are produced under one or more US patent. See www.abyrx.com/patent for more details.

ABYRX products are made in U.S.A.