

MONTAGE®

SETTABLE BONE PUTTY

Contains Hydroxyapatite/ β -TCP and 5% by weight Vitamin E Acetate

PREPARATION & APPLICATION TECHNIQUE GUIDE

STEP 1 PREPARE THE BONE

- Before opening the MONTAGE packaging, carefully inspect and clean the bone surface on which MONTAGE is to be applied.
- MONTAGE adheres best to bone that has been debrided to remove all undesired tissue (loose bone, periosteum, soft tissue) and any other foreign materials.
- Before applying MONTAGE, be sure to irrigate and dry to ensure a clean bone surface.




STEP 2 MIX MONTAGE

- Open MONTAGE packaging immediately prior to use and keep the components dry during the mixing process.
- For best ease of mixing, use the components of MONTAGE at room temperature (21°C / 70°F).
- Using sterile technique and dry gloves, vigorously knead and mix MONTAGE for at least 45 seconds, until MONTAGE is a single, uniform light tan color and consistency.



STEP 3 APPLY MONTAGE

- Using sterile technique and dry gloves and/or instruments, apply MONTAGE to the bone immediately following mixing.
- Be sure to apply enough pressure and spread MONTAGE into the pores of the bone so that it can interdigitate.

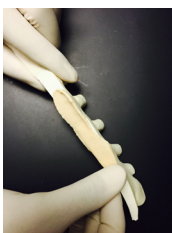


STEP 4 REMOVE EXCESS

- Remove any excess MONTAGE that is not necessary and/or impinging on adjacent soft tissues.
- If continued use of a mixed unit of MONTAGE is desired, it is very important to recombine any partitioned pieces of MONTAGE and/or apply MONTAGE immediately following mixing to ensure ease of application prior to hardening.
- Discard any unused putty.




IMPORTANT APPLICATION INFORMATION

- All bony surfaces on which MONTAGE is to be applied should be irrigated and free from undesired tissue (loose bone, periosteum, soft tissue) and any other foreign materials.
- All gloves and instruments that come into contact with MONTAGE should be dry; sterile technique should be followed from commencement of mixing to application to bone. The use of fresh gloves may be desirable.
- Using sterile technique and dry gloves, vigorously knead and mix MONTAGE for at least 45 seconds, until MONTAGE is a single, uniform light tan color and consistency.
- The ideal duration of application / working time (for adherence to bone and spreadability) is approximately within two (2) minutes of initiation of mixing.
- It is very important to remove from the bony site any excess MONTAGE that is not necessary and/or protruding into the adjacent soft tissues.

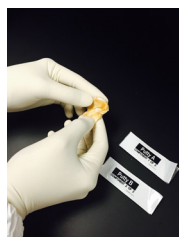

Combine and vigorously hand-mix the two putties for a minimum of 45 seconds

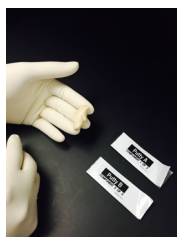

A fully-mixed Montage putty will have a uniform consistency and light tan color, and be minimally adherent to dry gloves and instruments



Instantly moldable



Instantly shapeable



Instantly cohesive



A fully-mixed Montage putty will harden within minutes



700 Fairfield Avenue, Suite 1 • Stamford, CT 06902 U.S.A. • 1.855.475.9175 • customerservice@abyrx.com • www.abyrx.com

© 2022 ABYRX, Inc.

ABYRX products are produced under one or more US patent. See www.abyrx.com/patent for more details.

ABYRX products are made in U.S.A.

D-000626 D