

FUZE/FIX

JONES SCREWS



A TRADEMARK OF
Fusion
ORTHOPEDICS

SCREW INFO

Self-Tapping / Self-Drilling

Screw allows easy insertion without pre-drilling.

Solid Screws

Durable titanium alloy material.

Torx Interface


The only interface that maximizes axial force and minimizes stripping is Torx.

Varying Thread Directions

Ability to maximize compression with counter-active thread forces.


Headless Compression Design

Allows the screws to be implanted with minimal risk of impingement or soft tissue irritation.

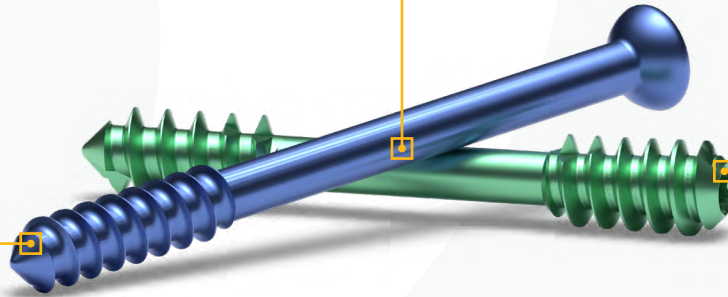

SIZING

	Diameter	Lengths
Blue	4.5mm	40mm - 70mm
Green	5.5mm	40mm - 70mm

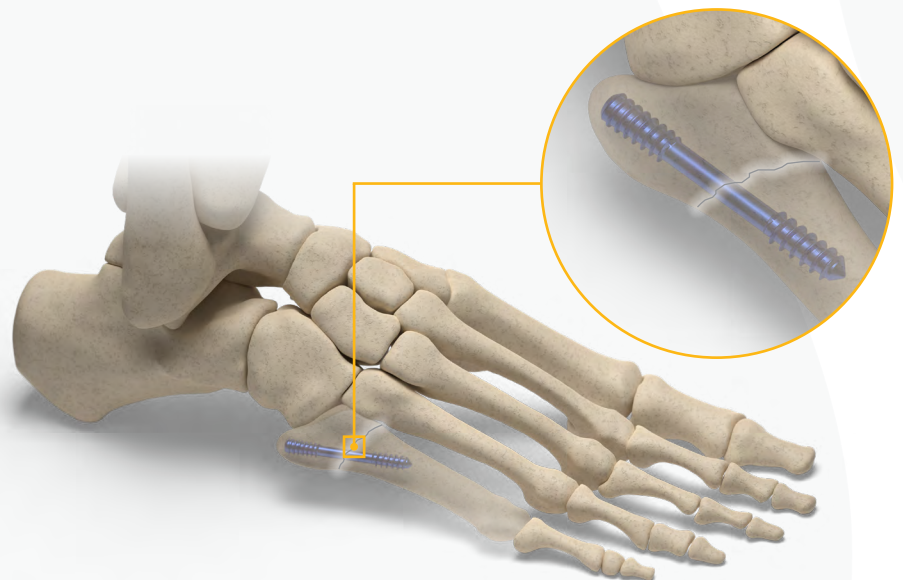

14 UNIQUE SCREW OPTIONS

FUSION ORTHOPEDICS USA, LLC.
 6859 E Rembrandt Ave #122
 Mesa, Arizona 85212 USA
 (800)-403-6876

Fusion
 ORTHOPEDICS

FDA Cleared.
 All Rights Reserved.
 FuzeFix is a Trademark
 of Fusion Orthopedics USA, LLC.
 SLS005 09132023

A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Surgeons should rely on their education and professional training, as the Fusion Instruments, implants, and techniques ("Products") augment traditional surgeries. Fusion Orthopedics, and its affiliates ("Fusion"), does not dispense medical advice and while Fusion recommends that surgeons be trained in the use of any particular Fusion Product before using it in surgery, surgeons must still rely on their professional training and expertise. The information presented is intended to demonstrate the breadth of Fusion's Products. A surgeon and patient must always refer to the current written surgical technique guide, fusionorthopedics.com and associated websites, package insert, product label, and instructions for use before considering and using any Fusion Product. Products and associated accessories may not be available in all markets because of Product availability, and Products are subject to the regulatory law and/or medical practices in individual markets. Please contact your representative if you have questions about the availability of Products in your area. Patients must solely rely on advice from their surgeon regarding the risks and all steps of a particular surgery with the utilization of Fusion's Products. As with any surgery, individual results and recovery time vary. Recovery time may range from months to years. Known risks associated with surgery are, but not limited to, injury, significant scarring, nerve damage, mal-union, long-term pain, implant pain, non-union, misalignment, loss of motion, bone shortening, infections, swelling, stiffness, bone fracture, and subsequent required amputation. These works, and its contents, are protected by U.S. and International copyright laws. Reproduction and/or distribution of these works, and its contents, without the express written consent of Fusion Orthopedics USA, LLC is prohibited. © 2023 Fusion Orthopedics USA, LLC